


VibraPEP®

Oscillatory PEP Therapy System (OPEP)

CombinedPEP® for atelectasis


Position 1, 2, or 3 creates CombinedPEP® characterized as continuous positive pressure above baseline with applied pressure and flow changes. These settings maintain appropriate positive expiratory pressure (PEP) that promotes collateral ventilation.

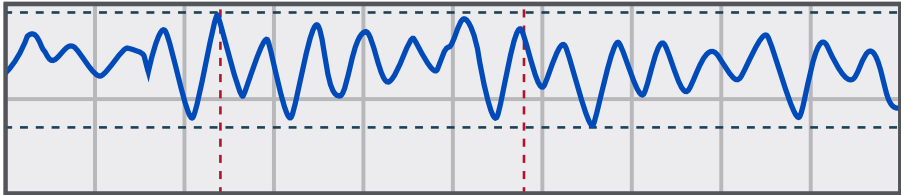

Oscillating Valve technology contorts based on the five therapy settings, allowing clinicians to tailor resistance to individual patient needs.

DynamicPEP® for bronchial secretions


Position 4 and 5 creates DynamicPEP®, erratic flow and pressure oscillations characterized by a slow increase from zero to maximum with an abrupt return to zero, assisting in the detachment and removal of bronchial secretions.

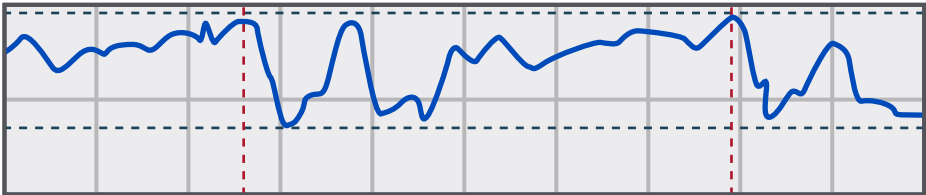



Placement of Filter Devices with the VibraPEP®

VibraPEP® Distal End

Distal end of the VibraPEP® allows an ISO 22mm OD filter to be affixed.

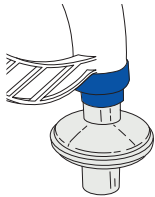

Adding a distal end filter can provide protection without compromising therapy.

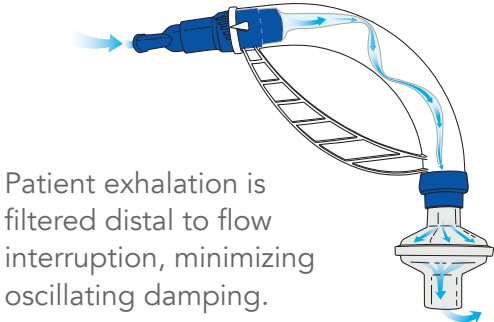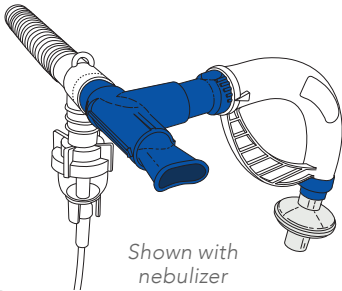

Flow Valve

Oscillating Valve Technology uses the patient’s expiratory flow to generate positive pressure and oscillations and is designed to function in any orientation



Removable Mouthpiece

Allows for adaptation to ISO connections

Comfort Handle

Integrated handle is easy to grip

Distal Filter Connection

The distal 22 mm connection permits attachment of a compatible respiratory filter

Valved Tee Adapter
Valved Tee Adapter with low resistance valve and ISO 22 mm OD x 18 mm ID Connection

Filter Connection
Filter attaches to the distal end of the VibraPEP®

Flow Valve
Provides for easy insertion, no introducer required


Therapy Selector
Clicks into place providing five treatment positions

Quick Start Instructions for Use:

- 1 Visually inspect the device for damage
- 2 Insert the mouthpiece
- 3 Inhale deeply through your nose, then hold your breath for 2 to 3 seconds
- 4 Exhale through the VibraPEP® for 3 to 4 times longer than your inhalation
- 5 Repeat this process for 15 breaths or as instructed by your healthcare provider
- 6 Perform 2 to 3 “huff coughs” to clear the airway

Instructions For Use with Nebulizer:

- A Assemble the nebulizer as described in the nebulizer supplied Instructions for Use
- B Connect the Valved Tee Adapter* and mouthpiece
- C Assemble Nebulizer to the 22 mm connector of the VibraPEP®
- D Turn on the nebulizer and inhale the medication through the mouthpiece and exhale through the VibraPEP® as in step 4

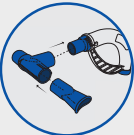

* The Valved Tee Adapter is an optional accessory Part #12218




with Nebulizer

Quick Start Cleaning Instructions:

Disassembly

- 1 Remove the Mouthpiece and End Cap
- 2 Rotate past position 1
- 3 Remove Therapy Selector and Valve from Tube
- 4 Insert Cleaning Aid into Valve during rinse Soap & Water

Cleaning Options

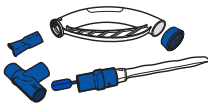

Soap & Water



Boiling Water
(10 minutes)



Dishwasher
(Top Shelf)



Air Dry
(Overnight)

Reassemble

- 1
- 2
- 3
- 4