

Easy Aerosol Delivery

For Your Tracheostomy Patients



Setting the standard

At Monaghan Medical we define the standard for special-need tracheostomy patients with the **AeroTrach Plus**® anti-static Valved Holding Chamber (aVHC). A simple and cost-effective system to quickly deliver medicine where it's needed most - straight to your patient's lungs.

Maximize Performance

The **AeroTrach Plus**® aVHC anti-static polymer chamber minimizes static charge and improves aerosol availability. The see-through chamber allows caregivers to visualize the medication dose in the chamber.

Quick Relief

The **AeroTrach Plus**® aVHC aerosol delivery is quick and easy – two important factors for any tracheostomy patient. Concentrations are delivered on inspiration only and are consistent – a benefit of the exceptional particle suspension time of the **AeroTrach Plus**® aVHC. As an additional benefit of the high particle concentrations, delivery is also cost-effective.

Expert Innovation

Monaghan Medical is renowned for its extensive research and development of new technology for aerosol delivery. We are proud to continue the tradition of innovation with the **AeroTrach Plus**® aVHC, an improvement in asthma relief for tracheostomy patients.



AeroTrach
Plus
anti-static Valved Holding Chamber



Of course it's from
monaghan™

AeroTrach Plus® anti-static Valved Holding Chamber (aVHC)

More efficient and comfortable medication delivery for your tracheostomy patients.


- **Comfort Right Fit™ Adapter**

Typical 15 mm adapters are designed to connect firmly in place. The **Comfort Right Fit™** adapter means the connection will fit tightly enough for an exceptional aerosol delivery, without excessive torque when connecting or disconnecting the device.



Valve and baffle system maximizes the delivery of fine particles emitted from standard pressurized MDIs

Unique, clear, Comfort Right Fit™ adapter ensures visibility

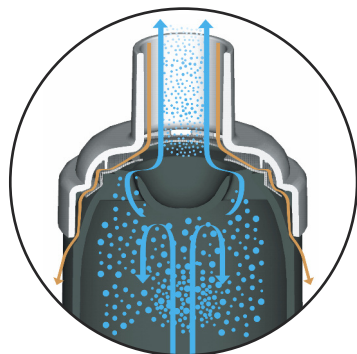

Flow Dynamic Design:
Improves respirable dose of medication.
Allows exhalation away from the patient.

- **Anti-Static Polymer**

Our anti-static polymer chamber minimizes static charge which improves aerosol availability.

- **Inhalation Valve**

With a protective design for low resistance and superior responsiveness, even the smallest inspiratory effort will be rewarded. In addition, the valve also features increased durability.

- **Exhalation Valve**

With our improved, more robust design, your patients will experience a low resistance to exhaled flow. Additional protection from valve tampering or damage has been built in.

- **MDI Adapter**

We've selected a polymer material for the MDI adapter. With this reinforced material, we've increased product durability.

Anti-static polymer material minimizes static charge which improves aerosol availability



Aerodynamic Chamber



For more information or a demonstration, contact your Monaghan Medical Sales Representative.

Ordering Information

Part #	Description	
52510	AeroTrach Plus® anti-static Valved Holding Chamber (aVHC)	cs/10



monaghan™

Monaghan Medical Corporation

5 Latour Ave., Suite 1600 • Plattsburgh, NY 12901 • **Customer Service 800-833-9653** • www.monaghanmed.com

™ and ® are trademarks and registered trademarks of Monaghan Medical Corporation or an affiliate of Monaghan Medical Corporation
© 2003, 2008 Monaghan Medical Corporation