




STRIVE® DUAL ZONE PEAK FLOW METER

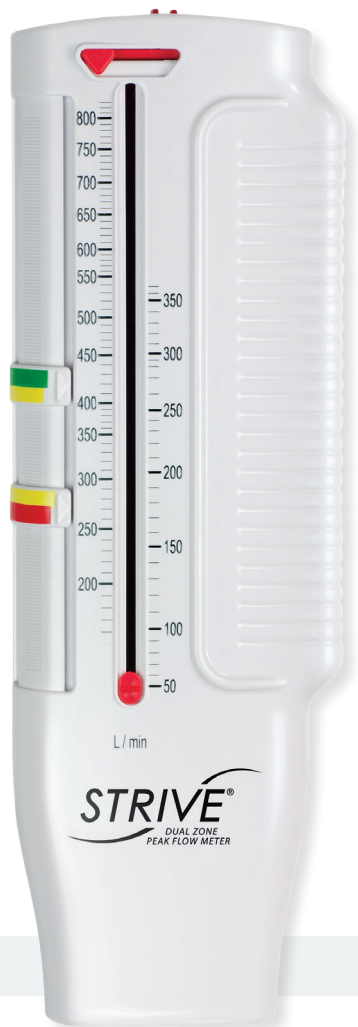

The Strive® Dual Zone Peak Flow Meter (PFM) is designed with both the clinician and patient in mind. For patients, the small, slim device is designed to fit in your pocket, ready to go anywhere. For clinicians, the Strive® PFM is a universal device with the versatility to provide accurate measurements for any patient population.

The Strive® PFM provides accurate measurements for patient populations from 5 years old through adult. Device can accommodate measurements from 50–800 L/min with the ability to measure two distinct zones: 50–350 L/min and 200–800 L/min.

Designed for individual or multi-patient use, the Strive® PFM is easy to clean for high-level disinfection and has optional disposable mouthpieces with one-way valves for clinical applications.

The Monaghan Advantage

- ❖ Portable device makes it easy to track day-to-day changes in peak expiratory flow
- ❖ Provides accurate measurements for patient populations from 5 years old to adult
- ❖ Measurements from 50–800 L/min accommodated
- ❖ Measures two distinct zones: 50–350 L/min and 200–800 L/min
- ❖ Easy to clean, top shelf dishwasher safe
- ❖ Not made or manufactured with natural rubber latex, phthalates, lead or PVC
- ❖ Designed and manufactured in North America

MEETS U.S. AND
INTERNATIONAL
TECHNICAL
STANDARDS

THE FACTS

The Strive® PFM - the one device that can simplify inventory needs and reduce costs.

Meets the most recent Technical Standards of Spirometry for portable peak flow meters established by the American Thoracic Society.¹

Meets International Standard for peak expiratory intended for the proper assessment of pulmonary function.²


REFERENCES

1. ATS/ERS Task Force: "Standardization of Lung Function Testing," 2005.
2. International Standard ISO 23747:2007: "Anaesthetic and Respiratory Equipment - Peak Expiratory Flow Metrics for Assessment of Pulmonary Function in Spontaneously Breathing Humans."



518.561.7330 phone
800.833.9653 toll free
518.561.5088 fax

5 Latour Avenue, Suite 1600
Plattsburgh, New York 12901
monaghanmed.com

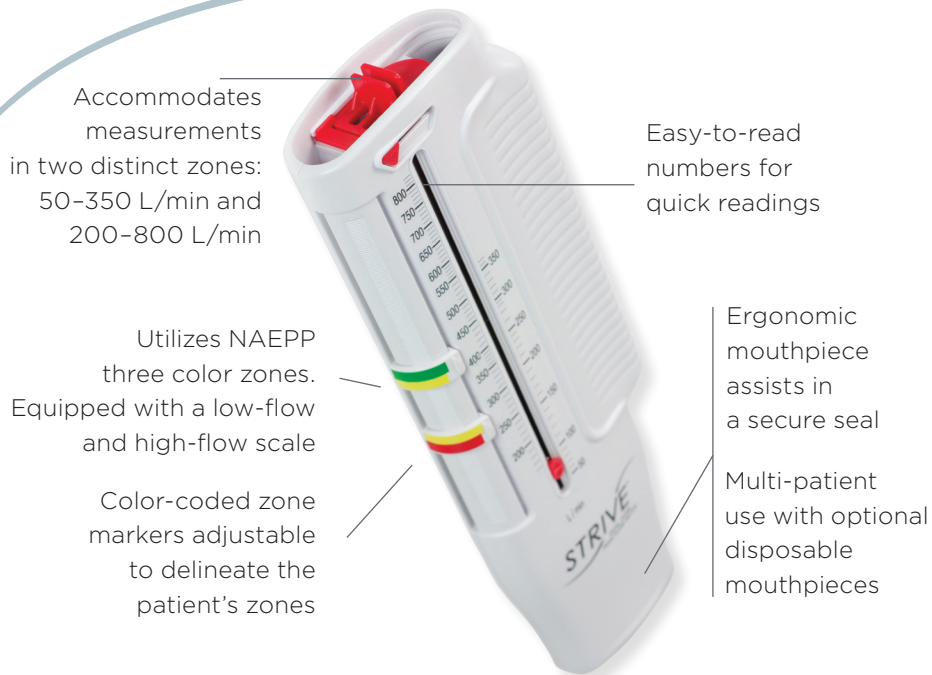

Accommodates measurements in two distinct zones: 50-350 L/min and 200-800 L/min

Easy-to-read numbers for quick readings

Utilizes NAEP three color zones. Equipped with a low-flow and high-flow scale

Ergonomic mouthpiece assists in a secure seal

Color-coded zone markers adjustable to delineate the patient's zones

Multi-patient use with optional disposable mouthpieces

The Strive® Dual Zone Peak Flow Meter meets the educational guidelines recommended by the National Asthma Educational Prevention Program.

The Strive® PFM combines a low flow scale and standard scale in one device, providing the option to replace two devices with one device, reducing the number of items to inventory. With only one device to choose from and one device to inventory, the Strive® PFM saves both time and money.

CLEANING IS EASY

- Designed for one year of use in the home, clinic or hospital environment
- Approved for the top shelf dishwasher cycle and high level disinfectant such as Cidex†

ORDERING INFORMATION

PART #	NAME/DESCRIPTION	UNITS PER CASE
96610	STRIVE® Dual Zone Peak Flow Meter (PFM)	10
90200	Disposable mouthpiece	200

Contact your Monaghan Sales Territory Manager for more information or a demonstration.