

# System



# **Choice of Power Source**

The EasyGoVac models give you the option of a Lithium battery for faster re-charge times and no memory, a lead acid battery or an AC unit for lightweight non-portable use.

### High Flow Rates

Even with filters, all models have high flow rates of between 21 and 26 LPM, with vacuum levels of 50 - 550 mmHg.

## **Canister Options**


All units are available with a dedicated 1,200 cc re-useable, autoclaveable canister or a standard 800 cc canister.

# Built-in Battery Charger

All battery models have a built-in charger, eliminating the need to carry a heavy external charger which can also be lost. Everything is inside the unit, to charge it, all you need to do is plug it in.

# Easy to Use

All models have an easy to read "power source indicator", "battery test indicator" and a "fully charged indicator" to make operation by the caregiver and/or Made in patient easier. A convenient carry bag is included with all models.

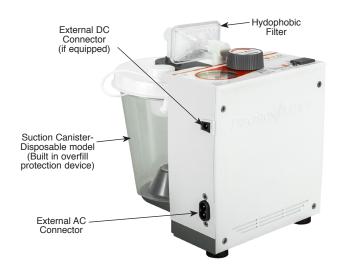



# **SPECIFICATIONS**

		PM66AC, PM66Li		PM66S	
		Disposable Bottle (800cc)	Reusable Bottle (1200cc)	Disposable Bottle (800cc)	Reusable Bottle (1200cc)
Dimensions:	Length	8.60 in. (21.8 cm)	7.45 in. (18.9 cm)	9.43 in. (23.9 cm)	8.28 in. (21.0 cm)
	Width	6.30 in. (16.0 cm)	6.30 in. (16.0 cm)	6.30 in. (16.0 cm)	6.30 in. (16.0 cm)
	Height	8.02 in. (20.4 cm)	7.38 in. (18.7 cm)	8.02 in. (20.4 cm)	7.38 in. (18.7 cm)
Weight:	(Device)			4.90 lbs (2.22 kg)	5.00 lbs (2.27 kg)
	PM66AC:	3.25 lbs (1.47 kg)	3.35 lbs (1.52 kg)		
	PM66Li:	3.70 lbs (1.67 kg)	3.80 lbs (1.72 kg)		
	(Shipping)			7.20 lbs (3.27 kg)	7.35 lbs (3.33 kg)
	PM66AC:	4.50 lbs (2.04 kg)	4.65 lbs (2.11 kg)		
	PM66Li:	6.00 lbs (2.72 kg)	6.15 lbs (2.79 kg)		

# Vacuum:

2 - 21 inHg (51 - 533 mmHg)

#### Flow

≥ 19 l/min @ 21.65 inHg (550 mmHg)

# **Battery Duration:**

2 hours

# **Battery Re-charge time:**

2 hours

## Flow:

26 LPM w/standard bacteria filter 21 LPM w/hydrophobic filter

# **AC power Source:**

100-240 VAC 5/60 Hz

#### Warranty:

2 years on unit

6 months on battery

Specifications are subject to change without prior notice.

# ORDERING INFORMATION

Part Number	Description	Part Number
PM66AC	EasyGoVac Aspirator (AC only)	502519-10
	includes Carry Bag	502954
PM66Li	EasyGoVac Aspirator with	503106
	Lithium Battery includes Carry Bag	507263
PM66S	EasyGoVac Aspirator with Sealed Lead Acid Battery includes Carry	506851
1 111003		507521
	Bag	507522
1955	Patient Suction Tubing	

Part Number	Description
502519-10	800cc Canister, (Case-10 ea.)
502954	Auto DC Power Cord
503106	AC Power Cord
507263	Carry Bag
506851	Inlet Filter (Hydrophobic)
507521	800cc Inlet Adapter Kit
507522	1200cc Canister

Precision Medical, Inc.

300 Held Drive • Northampton, PA 18067 USA

Phone: (+001) 610-262-6090 Toll-Free Phone: 1-800-272-7285 Fax: (+001) 610-262-6080 Toll-Free Fax: 1-800-352-1240

