


SALTER LABS®



Salter AIRE Plus® Compressor

Efficient aerosol delivery is an important part of managing the disease state of the patient. Efficient and convenient products mean delivering the aerosolized medication to the site in the shortest time while minimizing medication waste. Aerosol delivery with the Salter AIRE Plus® and Salter line of Nebulizers, including the new NebuTech® HDN®, helps to accomplish this and improve compliance.

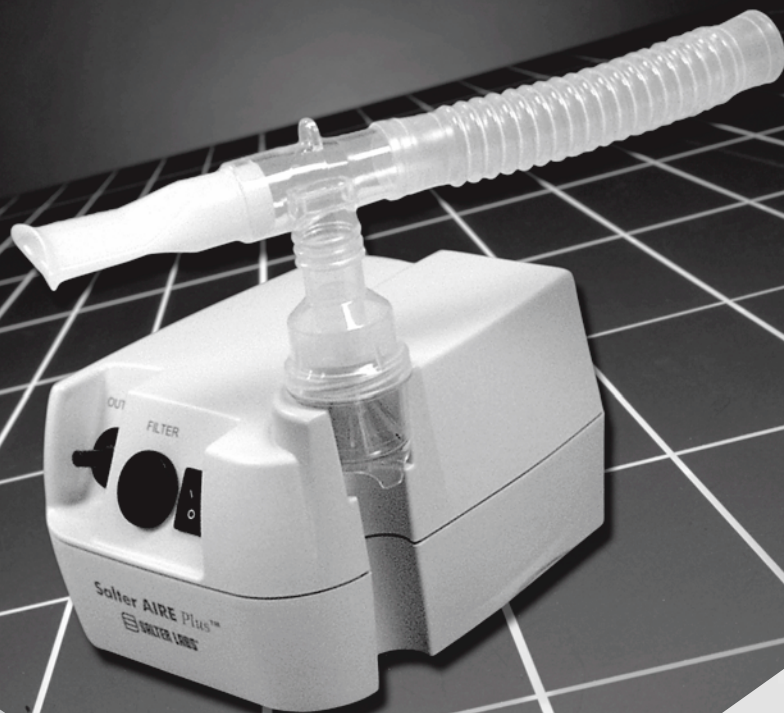

Features

- Compact
- Solid state components, 5 year warranty
- Lightweight
- Piston style compressor
- Quiet, low noise level
- Matched compressor and nebulizer
- Nebulizer docking port

Benefits

- Easy to use
- Durable and reliable
- Easy to carry, 4.0 lbs.
- Assures strong steady pressures and flows
- Provides unobtrusive treatment time
- 80% of particles in respirable range
- Added safety and ease of operation

Variety of uses. . . Hospital, Skilled Nursing Facility, Homecare, Clinic




Salter AIRE Plus® Compressor

A lightweight, portable, quiet compressor for powering aerosol generating medication nebulizers. It is engineered to provide long term, cost effective, rugged performance. The flow outlet fitting, fresh air inlet filter, and rocker style on/off switch are recessed to protect from accidental breakage. There is a convenient built in hand-grip for easy transport. A unique, tip proof nebulizer docking port is built into the case. An optional handy travel bag with extra storage for nebulizer kits is available.

Specifications

Salter AIRE Plus®, Model 8350 Compressor

Size:	6.7" L x 5.7" W x 4.4" H (17 cm X 14.5 cm x 11.2 cm)
Weight:	4.0 lbs. (1.8 kg)
Electrical requirements:	115 V, 60 Hz, 1.9 A
Compressor:	Oil free cast aluminum alloy integrated uni-piston
Free Flow	10 LPM, minimum
Liter Flow	6 LPM, minimum through a Salter Labs nebulizer
Pressure	25 PSI, minimum
Replaceable, external filter:	White, open cell foam
UL Classified to	UL 2601-1 / CAN/CSA C22.2 no.601.1 (See complete marking on product)

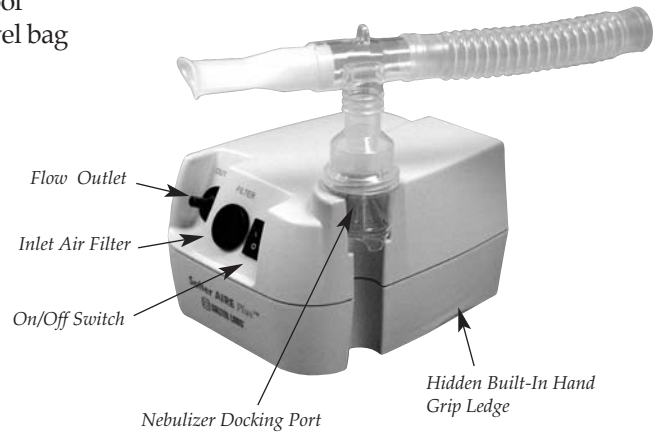

Pictured with 8900 series nebulizer

Ordering Information

	Part #
Compressor, Nebulizer with Tee Adapter, Mouthpiece and Reservoir Tube and 7' Tubing	8350-8900
Same as above, plus Compressor Carrying Bag	8350B-8900
Compressor, Nebulizer, Adult Elastic Strap Style, Aerosol Mask and 7' Tubing	8350-8924
Compressor, Nebulizer, Pediatric Elastic Strap Style, Aerosol Mask and 7' Tubing	8350-8906
Compressor, NebuTech® HDN® Nebulizer, (Reusable) and 7' Tubing	8350-8660
Compressor, NebuTech® HDN® Nebulizer, (Disposable) and 7' Tubing	8350-8960
Compressor, NebuTech® HDN® Nebulizer, (Disposable) Adult Elastic Strap Style, Aerosol Mask and 7' Tubing	8350-8984
Compressor, NebuTech® HDN® Nebulizer, (Disposable) Pediatric Elastic Strap Style, Aerosol Mask and 7' Tubing	8350-8966

Accessories	Part #
Salter AIRE Plus® Compressor carrying bag	8258
Replacement air filters, 5 each per package	8501
Replacement air filter cap, 1 each	8504



100 W. Sycamore Road - Arvin, CA 93203
 Phone: 661-854-3166
 1-800-421-0024 (CA only)
 1-800-235-4203 (U.S. & Canada)
 Fax: 661-854-3850
 Toll Free Fax: 1-800-628-4690 (U.S. & Canada)
 International Fax: 661-854-3850
 E-mail: salterlabs@us.salterlabs.com
 www.salterlabs.com

CE
0482

ISO 9001/EN 46001 & ISO 13485